


# High Frequency Wand

User's Manual

## LUCSUER®

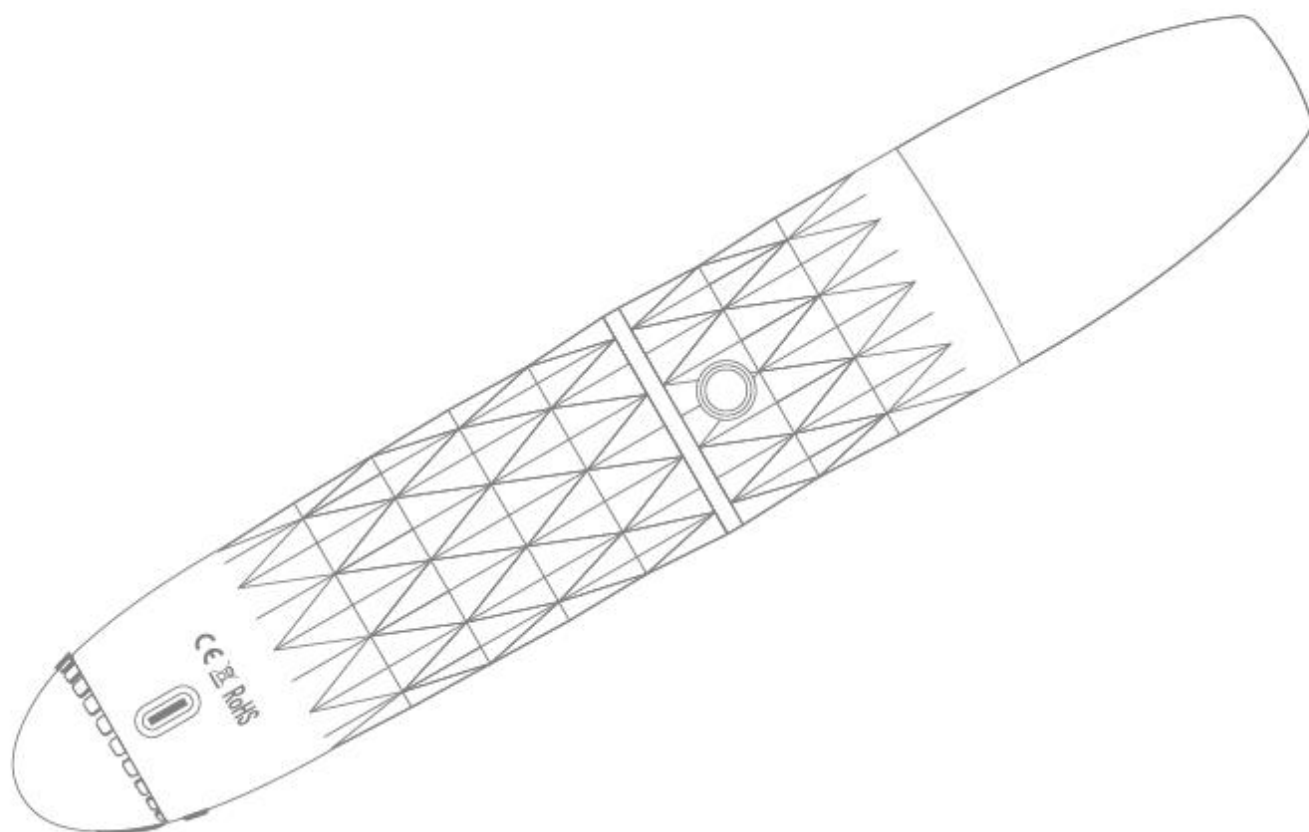

Thank you for choosing Our OxyGlow.  
Please read this manual carefully before use.  
Keep this manual for future reference.  
Product appearance may vary slightly from images.

# LUCSUER®

**Look Beautiful, Feel Luckier**

In today's fast-paced world, the pursuit of beauty and personal care is more than just a routine; it's a statement of self-confidence and well-being. At Lucsuer®, we understand the transformative power of looking and feeling your best. That's why we are proud to introduce our slogan:

**"Look Beautiful, Feel Luckier."**

"Look Beautiful, Feel Luckier" captures the essence of what Lucsuer® stands for. It's not just about outer beauty, but also the inner glow that comes from feeling confident in your own skin. Our products are designed to enhance your natural beauty, giving you that radiant look that turns heads and boosts your self-esteem.

At Lucsuer®, quality is paramount. We combine cutting-edge science with the finest ingredients to create products that are both effective and gentle. Our commitment to excellence means that when you choose Lucsuer®, you're choosing products that deliver results you can see and feel.

We invite you to experience the Lucsuer® difference. Discover how our products can help you look beautiful and feel luckier every day. Join our community of empowered individuals who are embracing their beauty and feeling more confident and fortunate in their lives.

🌐 <https://lucsuer.com>    📧 [amazon@lucsuer.com](mailto:amazon@lucsuer.com)    📧 [support@lucsuer.com](mailto:support@lucsuer.com)

# Contents

---

1. WARNINGS AND PRECAUTIONS(P1-3)
2. PRODUCT OVERVIEW(P4)
3. WHAT IS HIGH FREQUENCY  
AND HOW DOES IT WORK?(P5-6)
4. INSTRUCTIONS FOR INSERTING AND REMOVING  
GLASS TUBES?(P7-8)
5. HOW TO USE THE OXYGLOW(P9-12)
6. HOW TO MAINTAIN THE OXYGLOW(P13)
7. MULTI-PURPOSED ATTACHMENTS(P14)
8. BATTERY USE AND CHARGING(P15)
9. TECHNICAL SPECIFICATIONS(P16)
10. TROUBLESHOOTING GUIDE(P17)
11. FREQUENTLY ASKED QUESTIONS(P18-19)
- 12. WARRANTY AND CUSTOMER SUPPORT(P20)**

# **WARNINGS AND PRECAUTIONS**

---

## **General Safety**

- Do not use the OxyGlow on infants or children;
- keep the OxyGlow and attachments out of reach of children and pets;
- For those with diminished physical, sensory, or cognitive abilities, OxyGlow use requires supervision by an adult or caregiver;
- Do not use if you are pregnant or if you have a heart condition;
- Consult your doctor before using OxyGlow with any medical issue or seizure risk.

## **Medical Considerations**

- Avoid using OxyGlow with implanted devices like defibrillators, neurostimulators, pacemakers, or ECG monitors, or if you have metal implants such as braces or permanent retainers;
- Before use, consult your doctor if you have any medical device or implant;
- Do not use OxyGlow if you have cancer or suspect skin lesions; seek medical advice;
- Stop use and consult your doctor if you experience pain or discomfort.

# WARNINGS AND PRECAUTIONS

---

## Device Care and Maintenance

- Use only the official OxyGlow Tubes and parts provided with your purchase;
- Do not use the OxyGlow if any parts, including attachments, are visibly damaged. Damaged components may affect safety and effectiveness;
- Store the OxyGlow in a cool, dry location, avoiding extreme temperatures;
- Do not attempt to disassemble or repair the OxyGlow yourself. If the OxyGlow malfunctions, has been damaged, or exposed to water, cease use immediately and contact us.

## Usage Instructions

- Remove all metal objects, including jewelry, before use;
- Ensure the Tube is inserted before powering on the OxyGlow;
- Use Bare Hands. Do not wear any gloves;
- Always hold the Power Ring for better experience and performance;
- Avoid use on severely dry or damaged skin;
- Can be used on the head but not over the eyes;
- Limit use to 15 seconds per area and no more than 5 minutes total;
- For external use only; not intended for internal application.

# **WARNINGS AND PRECAUTIONS**

---

## **Electrical Safety**


- Avoid using the OxyGlow near flammable substances, including liquids and gases, and do not use it with oil- or alcohol-based products;
- Do not expose any part of the OxyGlow to water; keep it away from sinks, faucets, and other water sources to prevent damage;
- Keep the OxyGlow away from other electronic devices to avoid interference or damage.

## **Charging and Battery Safety**

- Use only the provided USB-C cable for charging, and adhere to the manufacturer's guidelines;
- Do not operate the OxyGlow while it is connected to the charger;
- If the charging cable or adapter is damaged, contact the manufacturer; do not attempt repairs.
- When the OxyGlow reaches the end of its life, dispose of the built-in battery in accordance with local regulations.

# PRODUCT OVERVIEW

---

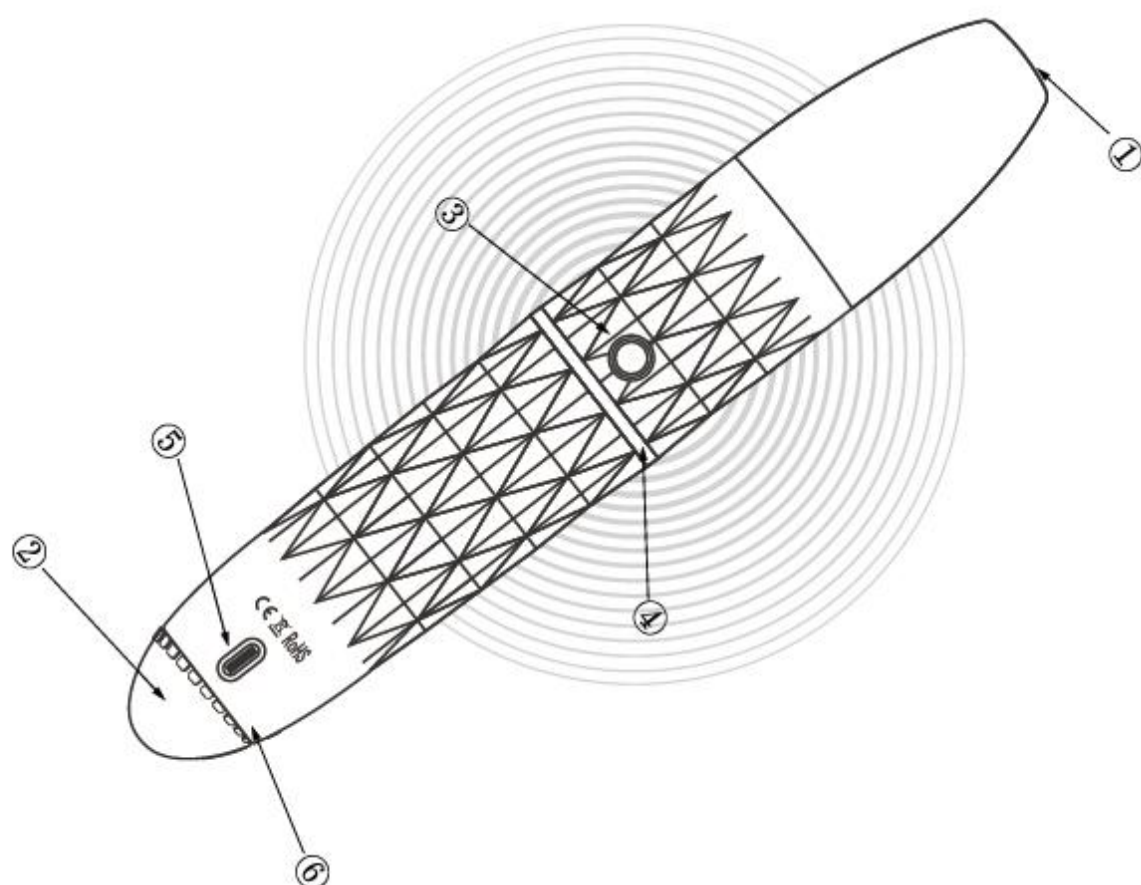

- ① **Electrode Chamber:** Generates high-frequency currents. Caution: Do not touch when on;
- ② **Intensity Dial:** Turns OxyGlow on and adjusts stimulation intensity;
- ③ **Safety Button:** Depress to activate current flow. Release to instantly cut power;
- ④ **Power Ring:** Hold for an ideal experience and performance when activated;
- ⑤ **USB-C Charging Port:** Enables convenient charging with any USB-C cable;
- ⑥ **Charging Indicator:** The indicator light is solid red while charging. It turns to solid green when the battery is fully charged.

## WHAT IS HIGH FREQUENCY AND HOW DOES IT WORK

---

### Understanding High Frequency

High Frequency is an advanced skincare technology that uses alternating electrical currents to treat various skin conditions. This method has been used in skincare for over a century.

### How High Frequency Works

1. **Device Composition:** Our OxyGlow Wand consists of a hand-held glass tube filled with xenon gas;
2. **Operational Principle:** Activation generates a mild electrical current from the wand;
3. **Skin Interaction:** The movement of the glass tubes across the skin produces gentle heat, mild electrical stimulation, and beneficial light energy;
4. **Ozone Generation:** The treatment also produces a minor amount of ozone, which boosts the efficacy of the skincare regimen.

### The Role of Xenon

We incorporate xenon in our OxyGlow Wand for its multifaceted benefits: it efficiently transforms electrical energy into light and heat; being inert, it ensures skin safety and prevents reactions with electrode materials; its stability extends the device's lifespan; and it emits specific wavelengths that effectively penetrate the skin.

## **WHAT IS HIGH FREQUENCY AND HOW DOES IT WORK**

### **Benefits of High Frequency Treatment**

The OxyGlow Wand delivers comprehensive benefits for both skin and hair health:

#### **1. Hair Loss Prevention and Scalp Health:**

It stimulates scalp blood circulation, enhances nutrient delivery to hair follicles, promotes denser hair growth, improves hair care product absorption, and regulates dandruff and oil production;

#### **2. Skin Brightening and Rejuvenation:**

It boosts microcirculation for a natural glow, accelerates cell metabolism, enhances the effect of brightening products, refines skin texture, and counters dullness from stress and environmental factors;

#### **3. Anti-Aging and Wrinkle Reduction:**

It encourages collagen and elastin production, reduces the appearance of fine lines and wrinkles, increases skin elasticity, relaxes facial muscles to minimize expression lines, and fosters a youthful radiance;

**4. Additional Benefits:** It provides acne treatment by eliminating bacteria and inflammation, improves lymphatic drainage to reduce puffiness and define facial contours, and enhances the efficacy of your skincare products.

## **INSTRUCTIONS FOR INSERTING AND REMOVING GLASS TUBES**

---

- 1. Unbox and Prepare the OxyGlow:** Remove the OxyGlow handle and the appropriate glass tube from the packaging. Always grip the glass tubes 3–5 cm above the metal base when handling;
- 2. Inserting the Glass Tube:** With the OxyGlow powered off, gently push the metal base of the glass tube into ①. Use steady pressure; some tubes may require more force to seat properly;
- 3. Never Apply Excessive Force:** Avoid forcefully jamming the glass tubes, as this could result in breakage and injury. Do not manipulate the tubes by their fragile tips;
- 4. Removing the Glass Tube:** Ensure the OxyGlow is powered off before removal. Grasp the glass tube close to the handle and lift it out directly;
- 5. Proper Handling Guidelines:** Avoid inserting any objects other than the official OxyGlow glass tubes. If any part of the wand or other objects becomes stuck, contact us immediately; do not attempt to remove it yourself.

## INSTRUCTIONS FOR INSERTING AND REMOVING GLASS TUBES

### IMPORTANT – READ THIS

#### THE GLASS TUBES ARE FRAGILE.

Do not handle them by the tip – they will break.  
Only handle them by the base when inserting and removing them from the OxyGlow handle.



**WRONG WAY**



**CORRECT WAY**

Warning: Always ensure the OxyGlow is powered off before inserting or removing glass tubes.


Always handle the glass tubes by holding them just above the metal base.

NEVER hold the glass tubes by the tips.

## HOW TO USE THE OXYGLOW


---

1. Begin by thoroughly cleansing your skin and patting it dry. Although high-frequency treatment can be applied to bare skin, using a water-based serum like High Energy Hyaluronic Serum can facilitate smoother wand movement for optimal results.

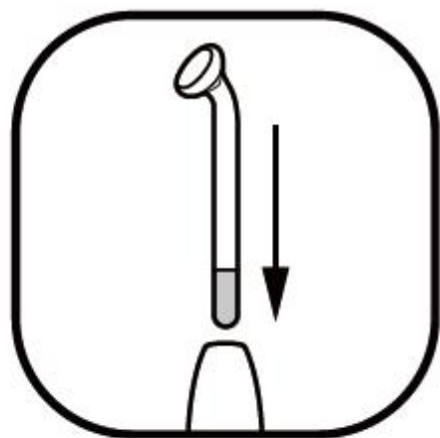

2. Select a glass tube and insert into ① securely.


3. Rotate ② clockwise to turn on. Adjust the intensity to your desired level. If you're a beginner, start with the lowest setting.

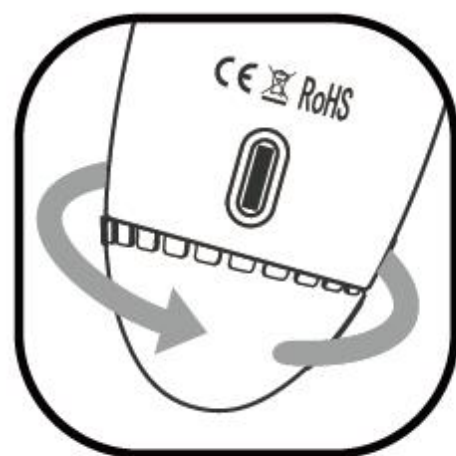

4. Hold the handle of your OxyGlow and place one finger at the base of the glass tube. This will prevent any small harmless static shocks when touching the tube to your skin.

## HOW TO USE THE OXYGLOW

---

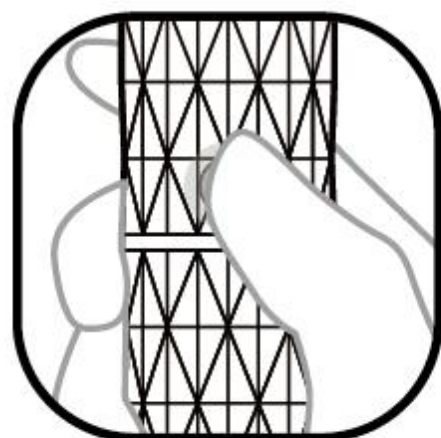

5. Press ③ to begin your session, and give a thumbs up to stop the OxyGlow when desired.



6. Test the intensity on your finger first. You should feel a gentle tingling, not a painful shock.

7. Upon activating the OxyGlow, you may hear a subtle buzzing or humming noise that increases with the power level. This sound is less noticeable when the glass tube is not in contact with your skin. Additionally, the tube may or may not emit a visible glow when in use, which is perfectly normal.

8. Apply 2-3 pumps of serum to one section of your face, working in small sections to avoid the serum from drying out.

9. With your finger on the glass tube, gently touch the tip to your skin, then remove your finger and glide the tube across your skin.

## HOW TO USE THE OXYGLOW

---



10. Use small, circular motions on your skin, spending no more than 15 seconds per area. Keep the tube moving; avoid prolonged contact on one spot. For targeted blemish treatment, tap the spot tube several times or hold it for 4-5s.

- Avoid the direct eye area, eyelids, lips, and sensitive areas.

11. After completing the treatment, massage the remaining serum into your skin, follow with your preferred moisturizer, or rinse off the serum after turning off the OxyGlow.



12. To turn off the OxyGlow after use, rotate ② counterclockwise, cleanse the OxyGlow, and store it properly.

**Caution:** For optimal device performance and gas tube longevity, restrict OxyGlow use to one session per day, not to exceed 5 minutes. This limitation helps maintain the gas integrity within the glass tubes.

## HOW TO USE THE OXYGLOW

---

"Begin with the lowest intensity setting for your initial use sessions. Gradually increase the level with each treatment as you become accustomed to the sensation, up to a point that is comfortable for your skin type and sensitivity. Remember, the system is effective even at the lowest settings, so there's no need to set the intensity too high to achieve optimal results."

## HOW TO MAINTAIN THE OXYGLOW

---

1. After your treatment, switch off the OxyGlow.
2. Clean the glass tubes with rubbing alcohol if necessary, and allow them to air dry completely.
3. Place the glass tubes back into their foam inserts. Ensure the tubes are completely dry before storing or using them again.

## MULTI-PURPOSED ATTACHMENTS



### **Mushroom(wide dispersal)**

Treats cheeks, forehead and jawline. Use this in conjunction with your pre-applied serums and creams to increase their absorption into the skin.

### **Comb(slight prickle dispersal)**

Stimulates hair follicles and promotes blood circulation. Diminishes hair loss while reducing dryness and dandruff.



### **90 Degree(intense and precise)**


Designed for acne, other skin imperfections, and smaller areas like crevices around the nose, lips etc.

**Inserting Glass Tubes:** Gently push the metal base of the glass tube into ①. Some tubes may slide in easily, while others may require firmer pressure.

**Removing Glass Tubes:** Ensure the OxyGlow is powered off; Hold the glass tube as close to the handle as possible (not the tip); Steadily pull the tube straight out of ①.

# BATTERY USE AND CHARGING

## Battery Warning

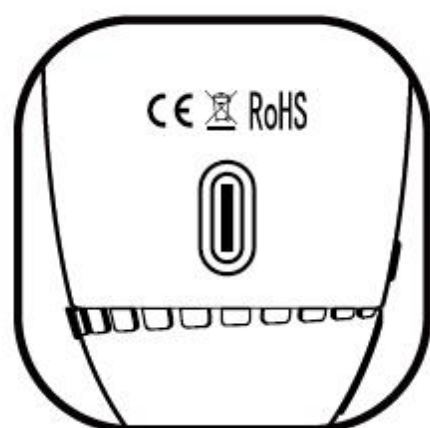

The OxyGlow Wand comes with a USB-C charging cable that works with any USB-C port. Find the charging port on the bottom of the wand's handle and connect the cable.

### Charging indicators:

- When charging: Solid red light.
- Fully charged: Light turns solid green

Fully charge the wand before first use. However, for the initial charge, avoid leaving it connected overnight to prevent overcharging.

The device can't work when charging. If you need to use it, please unplug the charging cable and use it.

Duration	2.5hours
Recharge	120minutes

# TECHNICAL SPECIFICATIONS

## Product Specifications

- Direct Current
- Power Mode: Battery
- Variable Frequency Dial Gauge
- Built-In Safety
- Materials:
  - PPV Plastic Hardware
  - Copper Tubes
  - Glass Rods
- 3 Glass Electrode Attachments.
  - 1 Comb Tube: 6 in.(L)x 1.25 in.(W)
  - 1 Mushroom Tube: 6 in.(L)x 1 in.(W)
  - 1 90 Degree Tube: 5.5 in.(L)x 1 in.(W)
- Voltage: 7.4V
- Current: input 1000mA output 0.2mA
- Input Power: 5W
- Dimensions: 13 in X 8 in X 2 1/3 in (33cm X 21 cm X 6cm)

### Frequency Range

30-1000MHz

### Net Weight

250g(0.55lbs)

### Gross Weight

550g(1.20lbs)

## Product Information

Product Name :	High frequency facial device
Model:	OG2-V1
Rated Voltage:	7.4 volts
Rated Current:	0.46A
FCC ID:	not applicable
Country of Origin:	Made in China
Manufacturer:	Dongguan ANE Crafts Co.,Ltd
Address:	14 Minghua Road, Shijie, Dongguan, Guangdong, P.R.C



# TROUBLESHOOTING GUIDE

Problem	Possible Cause and Solution
Device won't turn on	<ul style="list-style-type: none"><li>• Low battery: Charge the OxyGlow</li><li>• Power switch: Ensure it's turned on</li></ul>
Charging issues	<ul style="list-style-type: none"><li>• Check USB-C cable connection</li><li>• Try a different cable or power source</li></ul>
Device stops unexpectedly	<ul style="list-style-type: none"><li>• Low battery: Recharge</li></ul>
No glow or weak glow	<p>Please note that larger glass tubes like the mushroom tube may not produce any visible glow at all – this is normal and does not impact your results.</p> <p><b>1. Ensure close proximity:</b> The accessory is designed to emit light when in close contact with your skin.</p> <p><b>2. Adjust intensity:</b> Use ② to fine-tune the intensity to your preference. Turn clockwise to increase.</p> <p><b>3. Troubleshooting:</b> If the tube is damaged or the inert gas is depleted, the accessory will not function properly. Please visit <a href="https://lucsuer.com">Lucsuer.com</a> for assistance or to purchase a replacement.</p>

# FREQUENTLY ASKED QUESTIONS

---

## **Q1: How often and for how long should I use OxyGlow?**

Use OxyGlow no more than once every other day, limiting sessions to a maximum of 5 minutes. While treatment times may vary by skin type and concerns, exceeding these guidelines does not improve results and could cause redness. For optimal outcomes, follow the manual's recommendations and moisturize after each use to maintain skin hydration.

## **Q2: How quickly will I see results?**

Results may vary. Some users notice improvements such as better skin tone and reduced acne after a single session, while others may take several treatments. Since results are cumulative and dependent on individual skin types and concerns, consistency is key.

## **Q3: Should I use a conductive gel or serum on my skin during treatment?**

Apply a water-based serum to enhance the treatment. Avoid oil- or alcohol-based products as they are not suitable for use with OxyGlow.

## **Q4: Should I avoid any ingredients like retinol or glycolic acids?**

Avoid using exfoliating acids or irritating ingredients like retinol on days you plan to use OxyGlow. Schedule such treatments for alternate days to prevent skin irritation.

## FREQUENTLY ASKED QUESTIONS

---

### **Q5: Should I feel the electrical current on my skin?**

A mild tingling sensation is normal while using OxyGlow. Adjust the intensity using ② to match your comfort level. The OxyGlow remains effective at low settings, eliminating the need for high power. Consistent use as per guidelines is crucial. If you experience a 'shock' sensation, refer to the manual for correct OxyGlow application and finger placement.

### **Q6: Does high frequency kill bacteria?**

Yes. High-frequency treatments oxygenate the skin, targeting *P. acnes* bacteria, enhancing circulation, aiding detoxification, and reducing redness. They also help in pore minimization, anti-aging, and managing persistent acne.

### **Q7: Can I use this device immediately after bathing?**

For optimal performance and safety, allow your skin to dry completely before use. This ensures the OxyGlow operates effectively and reduces potential risks.

### **Q8: Is there a warranty for this device?**

This OxyGlow comes with a one-year limited warranty. We will repair or replace it at no charge if manufacturing defects occur under normal use.

# CUSTOMER SUPPORT AND WARRANTY

---

## Warranty Information

Lucsuer® warrants this OxyGlow against defects in materials and workmanship under normal use for a period of one year from the date of retail purchase by the original end-user purchaser.

1. **What is covered:** Defects in materials and workmanship under normal use.
2. **What is not covered:** Damage caused by accident, abuse, misuse, or improper maintenance; Products with lot numbers that have been removed, altered, or defaced; Normal wear and tear.

Simply scan the QR code right to start your **FREE WARRANTY** journey. Our customer support team is ready to assist you with any questions you may have.



## Customer Support

If you have any questions or concerns about your OxyGlow, please reach out to our customer support team for assistance:

- **Email:** [support@lucsuer.com](mailto:support@lucsuer.com)
- **Website:** <https://lucsuer.com/>